



**Meet the Alpha**

# THE GAME CHANGER IN PATIENT CARE

Move beyond the limitations of BMI and gain medically validated insights into your patient's body composition in just seconds – as easy as stepping on a scale, as precise as body composition gold standards!

**Helping you help.**





seca mBCA Alpha

# MEDICAL INSIGHTS THAT CHANGE LIVES

Get a clinically validated assessment of your patient's body composition – in just seconds. With clear, easy-to-understand visuals and meaningful metrics like skeletal muscle mass, visceral fat, and fat mass, you gain the insights you need to guide therapy, tailor interventions, and make informed decisions far beyond BMI and guesswork.

[seca.com/mbca-alpha](https://seca.com/mbca-alpha)



## The Key to Proactive Patient Management


True prevention starts with body composition analysis. When integrated into routine check-ups, it reveals critical developments – such as skeletal muscle loss or increased fat mass – at an early stage. Take action early and help patients build a strong foundation for lifelong health.



seca analytics 125

## See What Others Miss

Measure body composition with the seca mBCA Alpha and access results easily via your browser. Comprehensive parameters for adults – and dedicated ones for children – are at your fingertips. Powered by seca analytics 125, our continuously updated, cloud-based software that keeps you at the forefront of bioimpedance analysis.

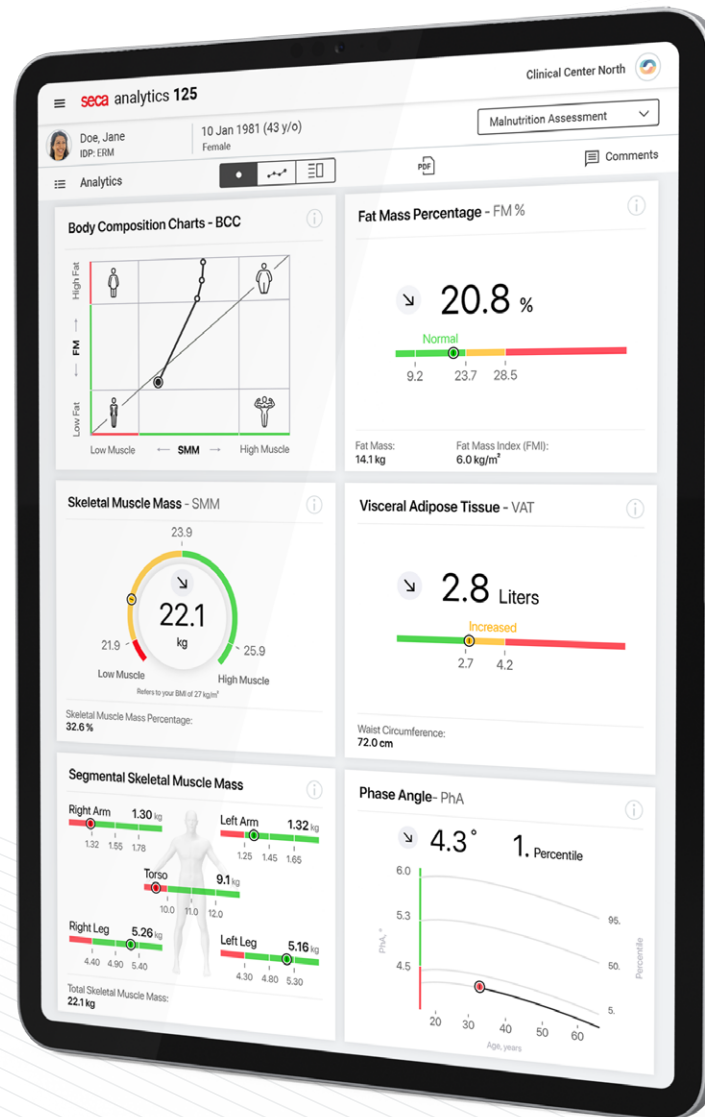

seca myAnalytics

## Put Results in Your Patients' Hands

With our app, your patients can access their body composition data anytime, anywhere. They can track progress, visualize goals, and stay engaged – right from their phone. Because motivation grows when progress is visible!



## Highest level of data protection



- Pioneer in data protection with ISO 27001 and ISO 27701 certification
- Effective protection of personal data according to GDPR
- System architecture developed from the ground up according to a security-by-design concept






## Validated Against Gold Standards


The seca mBCA Alpha sets a new benchmark in clinical body composition analysis. Its muscle mass readings are validated against whole-body MRI, and fat mass is benchmarked using the 4-compartment model, both gold standard in body composition analysis. Our validation process covers a wide range of patient groups – including children, adults over 65, individuals with a BMI over 30 kg/m<sup>2</sup>, and diverse ethnic backgrounds.

This comprehensive scientific foundation ensures your results are not just accurate, but clinically meaningful – backed by peer-reviewed research and ready to support your decisions.

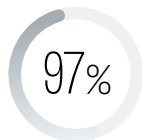

### Skeletal muscle mass

97% accuracy compared to whole-body MRI

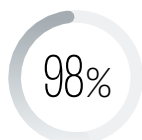

### Fat-free mass

98% accuracy compared to 4C model

## Book an appointment now

### Australia

☎ 1800 675 432  
✉ [enquiries@ecomед.com.au](mailto:enquiries@ecomед.com.au)

### New Zealand

☎ 0800 300 421  
✉ [sales@ecomед.co.nz](mailto:sales@ecomед.co.nz)



## seca mBCA Alpha

Medical device	Medical Device Regulation (MDR) certified
BIA measuring method	8-point Bioelectrical Impedance Analysis
BIA measuring frequency	1; 2; 5; 10; 20; 50; 100; 200; 500 kHz
BIA measurement segments	Right arm, left arm, right leg, left leg, right half of body, left half of body, torso
BIA measuring electricity	100 µA
BIA measuring period	24 seconds
Age	5 - 99 years
Body height	> 130 cm / > 51 inch
Software	seca analytics 125
Capacity	300 kg
Graduation	50 g
Dimensions (WxHxD)	600 × 1,233 × 653 mm
Platform (WxHxD)	600 × 70 × 467 mm
Net weight	20.8 kg
Display	4.3" tiltable and rotatable touchscreen
Power supply	Power adapter
Interfaces	Ethernet, USB for barcode scanner, Wi-Fi
Accessories	+ Barcode scanner mount seca 463 + Barcode scanner seca 473

Order code **ON IGM ANT NN**

- BIA platform seca mBCA 555
- BIA handle seca mBCA 545

For more information on our BIA solutions, visit us here:



[seca.com/mbca](https://seca.com/mbca)

**seca asia pacific**  
Lot B, 25-02, Q-Sentral  
Jalan Stesen Sentral 2  
50470 Kuala Lumpur  
Malaysia

☎ +60 173 711 973  
✉ [info.aspac@seca.com](mailto:info.aspac@seca.com)

**seca – Ecomed Trading pty ltd**  
5/30 Foundry rd, Seven Hills  
2147 NSW Australia

☎ 1800 675 432  
✉ [enquiries@ecomед.com.au](mailto:enquiries@ecomед.com.au)

**seca – Ecomed Trading pty ltd**  
PO Box 91044  
West Auckland  
1142 New Zealand

☎ 0800 300 421  
✉ [sales@ecomед.co.nz](mailto:sales@ecomед.co.nz)